Application and safety instructions for MEISINGER

Implant Removal Kit Implant Extraction Kit

CONTENT

Language	Page
English (EN)	2-8
French (FR)	9-16
German (DE)	17-23
Italian (IT).....	24-31
Polish (PL).....	32-39
Spanish (ES).....	40-47
Ukrainian (UA)	48-55

English (EN)

Application and safety instructions for MEISINGER

Implant Removal Kit Implant Extraction Kit

Item no.: BIR00, BIR01 and BIE00

The Implant Removal Kit enables simple and safe removal of implants inserted into the jaw. Specially developed threaded pins reliably connect the implant and the extraction instrument, allowing the torque to be transferred in a controlled manner and the implant to be extracted without damaging the tissue. Depending on the situation, a new implant with the same diameter can be inserted after the implant has been successfully removed.

The Implant Extraction Kit enables gentle, guided explantation of osseointegrated implants using precisely selected trephines. The simple guidance system allows the trephine diameter to be optimally adapted to the implant, which protects the surrounding bone tissue. The kit offers a variety of sleeves, threaded pins, and trephines for individual selection.

INSTRUCTIONS FOR USE

The Bone Management System Implant Removal Kit and the Bone Management System Implant Extraction Kit are indicated for the removal of dental implants from the jawbone that are blocked during implantation and/or whose surrounding bone has been damaged.

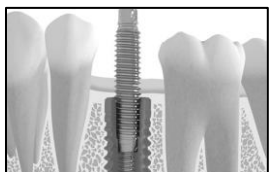
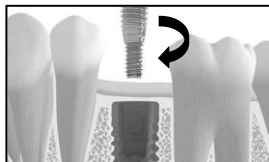
The Implant Removal Kit is indicated for the removal of implants where the degree of osseointegration allows removal by unscrewing within a risk-free torque range. In cases of advanced osseointegration of the implant, trepanation and use of the Implant Extraction Kit is indicated. The selection of the implant removal technique with the lower risk is made on a case-by-case basis by the user.

APPLICATION

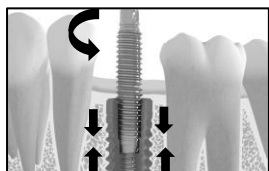
The maximum permissible torques for use and the recommended speeds are summarized in the article overview under "Contents." These must be observed to ensure proper and safe use.

Explantation of the implant with the Implant Removal Kit

First, select a threaded pin from the Implant Removal Kit that fits into the internal thread of the implant to be removed. The selected threaded pin can now be turned into the implant (clockwise) using the SWRA1 screwdriver.

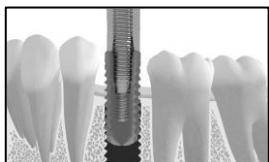
The ratchet is then attached to the screwdriver, allowing the threaded pin to be screwed into the implant with a torque of up to 30 Ncm.

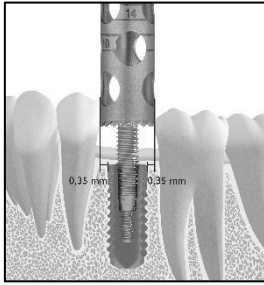
After removing the ratchet and screwdriver, a suitable extraction tool is selected. The ratchet adapter RAITR is placed on the extraction tool and then turned counterclockwise by hand as far as possible onto the threaded pin.

Using the ratchet, which is placed on the ratchet adapter, a torque of up to 100 Ncm can now be applied in the extraction direction (counterclockwise) until the implant rotates and is extracted.

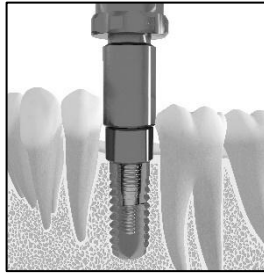
If the implant is difficult to remove, a rose bur can be used to remove a minimal amount of bone from the upper part of the implant. To do this, the extraction tool must be removed. The previous steps for explanting the implant must then be repeated.

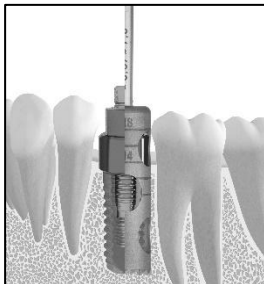
Explantation of the implant with the Implant Extraction Kit

First, select the appropriate trephine from the Implant Extraction Kit according to the diameter and length of the implant to be removed. Ensure that the inner diameter of the trephine is at least 0.7 mm larger than the widest part of the implant

Next, select the threaded pin that fits into the internal thread of the implant. Screw the corresponding sleeve, which is matched to the diameter of the selected trephine, onto the threaded pin. It is important to note that the sleeve should only be screwed onto the threaded pin to the extent that the SWRA1 screwdriver can be placed completely onto the threaded pin. The threaded pin is screwed into the implant together with the sleeve and tightened with the SWRA1 screwdriver (clockwise) to a maximum torque of 30 Ncm. The sleeve is then fixed in place using the SWRA2 screwdriver with a torque of 35–40 Ncm, ensuring that it rests correctly on the implant shoulder.

During extraction, the trephine is guided axially over the screwed-on sleeve. An optimum speed of 300–500 rpm must be maintained and sufficient cooling must be ensured. The depth of the trepanation is controlled by the user depending on the implant length and taking into account the adjacent anatomical structures. If the trephine hits the threaded pin, the procedure can either be continued with a longer trephine, or the threaded pin and sleeve can be removed in order to continue the trepanation to the desired depth. In this case, the trephine is guided through the previously created notch, while continuing to ensure that it is guided along the correct axis.

After explantation of the implant

To remove the threaded pin or the extraction instrument/sleeve from the explanted implant, the implant must be held in place with an auxiliary tool, such as pliers. The extraction tool/sleeve can now be turned out clockwise using the ratchet adapter/screwdriver SWRA2 and the ratchet.

After the extraction tool/sleeve has been removed, the SWRA1 screwdriver is attached to the threaded pin and the ratchet is attached to it. The threaded pin can then be unscrewed from the implant by turning it counterclockwise using the screwdriver and ratchet.

Sterilization

The instruments are supplied non-sterile and must be reprocessed (cleaning, disinfection, and sterilization) before use. For proper reprocessing of the items, please refer to the instructions for reprocessing (cleaning, disinfection, and sterilization) of medical devices from Hager & Meisinger GmbH at <https://www.meisinger.de/services/downloads>.

Disposal

When disposing of the instruments (after use or after the specified shelf life has expired), ensure that the product is disposed of in the waste for biological hazardous substances. All packaging components are disposed of in accordance with national regulations (e.g., dual waste system).

Intended use

The products are indicated for the removal of dental implants from the jawbone that are blocked during implantation and/or whose surrounding bone has been damaged.

Intended patient group

Taking into account the contraindications, the products can be used in patients of any gender.

Intended user profile

Professional users/specialist users in the field of dental implantology.

Contraindications

In general, general medical as well as local, absolute, and relative contraindications for dental surgical procedures must be observed. After successful explantation and when planning to implant a new implant, the contraindications for the implantation of an endosseous dental implant should also be observed. Information on these can be found in the instructions for use of the selected implant system from the respective implant manufacturer. There are no specific contraindications for the use of the instruments contained in this kit.

Care must be taken to protect the anatomical structures (safety distance of at least 2 mm) and to avoid damaging the adjacent teeth/tooth roots (risk of damage, infection/dehiscence).

ASPIRATION RISK: Care must be taken to ensure that the instruments are secured against aspiration during treatment.

When selecting the size of the threaded pin, the exact size of the internal thread of the implant to be removed must be known. Otherwise, incomplete fastening of the threaded pin in the implant may result in damage or breakage of the threaded pin before the implant can be removed. If the exact size is not known, it is recommended to start with the largest threaded pin and work down until the correct threaded pin size for the implant is found. When selecting the size of the tapping tool, it is recommended to choose a size that is larger than the diameter of the upper part of the implant. If a smaller size is used, the removal torque cannot be properly absorbed, making it difficult to remove the implant. When using the instruments, the maximum torques must not be exceeded. Otherwise, the threaded pin or removal tool may break or bend, the implant may splinter, or the bone may break. If a torque of 60 Ncm or more is applied during implant removal, the implant must also be cooled sufficiently to prevent overheating. **The inner diameter of the trephine must be at least 0.7 mm larger than the maximum diameter of the implant. The sleeve must then be selected to match the diameter of the selected trephine.**

The instruments may bend or break when reused or used multiple times. It is therefore the responsibility of the professional user to check the instruments for wear and tear after each use and to assess whether they can be reused after reprocessing.

The instruments in the Implant Removal Kit and **Implant Extraction Kit** were developed for implants from Meisinger Implants GmbH and myplant GmbH. The use of the instruments for other implant systems is the responsibility of the professional user.

The following applies to all milling and drilling procedures: In order to reduce the risk of bone overheating and thus necrosis, the respective instrument must be used intermittently with low pressure and under constant cooling with sterile physiological saline solution.

Serious incidents related to the product must be reported to the manufacturer and the competent authority of the Member State in which the user of the instruments for other implant systems is located.

General information

Please also note the general application and safety instructions for MEISINGER products in the medical field and the instructions for reprocessing (cleaning, disinfection, and sterilization) of medical devices from Hager & Meisinger GmbH at <https://www.meisinger.de/services/downloads>.

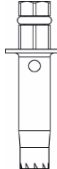
	REF	MD	LOT			CE/CE 0044
Manufacturer	Item number	Medical device	Batch number	Non-sterile	Caution	Conformity mark
			 ifu.meisinger.de	UDI	CH REP	UK RP
Do not use if packaging is damaged	Date of manufacture	Packaging unit	Consult the electronic instructions for use	Unique device identification	Swiss authorized representative	UK authorized representative

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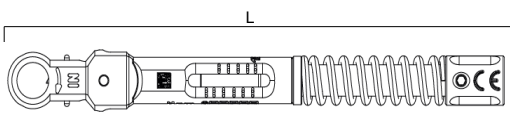
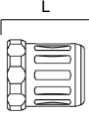
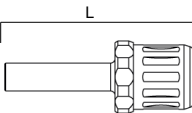
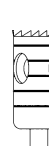
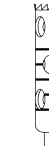
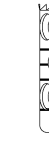
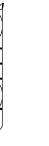
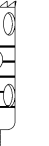
	Ratchet adapters	Screwdriver for threaded pins
		
	1	1
Fig.	RAITR	SWRA1
Length mm	8,5	20,5
		

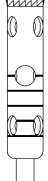
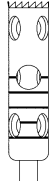
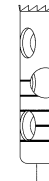
Contents BIE00

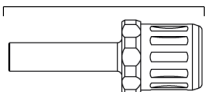
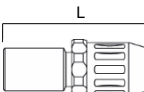
	Threaded pin				
					
	2	2	2	2	2
Fig.	GS6KT	GS6KT	GS6KT	GS6KT	GS6KT
Size	M1,6	M1,8	M2,0	M2,5	UNF1-72
Length mm	15,5	15,5	15,5	15,5	15,5
Max. torque Ncm	30	30	30	30	30
Color	Blue	Yellow	Orange	Red	Green
					

	Sleeve				
					
					
	1	1	1	1	1
Fig.	BIE40	BIE45	BIE50	BIE55	BIE60
Size	040	045	050	055	060
Length	6.0	6.0	6.0	6.0	6.0
Max. torque Ncm	40	40	40	40	40
					

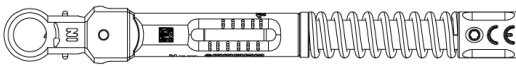
	Trephine						
							
	1	1	1	1	1	1	1
Fig.	229	229	229L	229L	229L	229L	229L
Shaft ¹	205	205	205	205	205	205	205
Size ²	055	060	040	045	050	055	060
Length mm	10,0	10,0	14,0	14,00	14,0	14,0	14,0
 mm	5,5	6,0	4,0	4,5	5,0	5,5	6,0
 mm	6,5	7,0	5,0	5,5	6,0	6,5	7,0
Optimal thread pitch rpm	300-500	300-500	300-500	300-500	300-500	300-500	300-500
Maximum thread pitch rpm	2.000	2.000	2.000	2.000	2.000	2.000	2.000
	 0044						

¹ 205 = RA L, ² Largest workpiece diameter in 1/10 mm  Inner diameter  Outer diameter

	Trephine			
				
	1	1	1	1
Fig.	229XL	229XL	229XL	229XL
Shaft ¹	205	205	205	205
Size ²	040	045	050	060
Length mm	18,0	18,0	18,0	18,0
 mm	4,0	4,5	5,0	6,0
 mm	5,0	5,5	6,0	7,0
Optimal thread pitch rpm	300-500	300-500	300-500	300-500
Maximum thread pitch rpm	2.000	2.000	2.000	2.000
	CE 0044			

	Screwdriver for threaded pins	Screwdriver for sleeves
		
	1	1
Fig.	SWRA1	SWRA2
Length mm	20,5	19,0
		

Optionally available:

	Torque ratchet	
		
	L	
	1	
Fig.	CARA3	
Length mm	109,0	
	CE	

¹ 205 = RA L, ² Largest workpiece diameter in 1/10 mm  Inner diameter  Outer diameter

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